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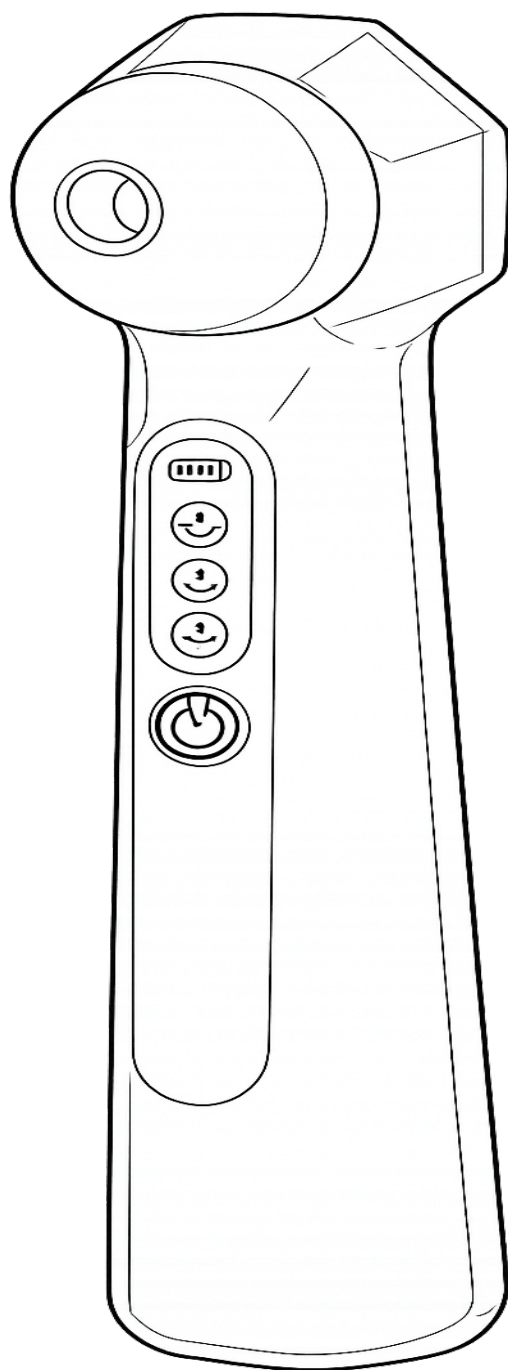
Dear Sir or Madam, thank you for purchasing our product!

Before using the product, please read carefully with the instructions below for proper use of the product.

Please keep this manual for future reference and following its recommendations, because failure to follow it may pose a threat to life or health.

Introduction to the user manual

This user manual has been prepared to help you familiarize yourself with the proper and safe use of the device. It contains important instructions on operation, maintenance, and safety rules, which, if followed, will ensure the smooth operation of the device and extend its service life.

Before first use, please read the instructions carefully. The information contained herein will also help you properly store the device, especially if it will not be used for an extended period.

Following the recommendations below will allow you to enjoy trouble-free and safe operation of the product.

ATTENTION!

Using the device in a manner inconsistent with the information contained in the manual may result in the risk of bodily injury.

Application and description of the device

The Blackhead Vacuum Cleaner is a modern, handy device designed for deep facial cleansing. Using vacuum technology, it effectively removes impurities trapped in pores, excess sebum, and blackheads. The device is equipped with interchangeable attachments and adjustable suction intensity, allowing it to be adapted to various skin types, including sensitive skin. Built-in LED lighting facilitates precise treatment.

This device is intended for home use only and is not intended for medical or commercial use, and is not intended to replace medical therapy. The device is intended for indoor use in dry rooms. It is not intended for medical, therapeutic, or commercial use.

For safety and CE certification reasons, this product must not be rebuilt or modified in any way. Using the product for purposes other than those described may result in damage to the product. Improper use may also result in hazards such as short circuits, fire, electric shock, etc.

WARNING!

- Danger of improper use!
- Using the device in a manner other than its intended purpose and/or using the device in any other way may involve various risks.

-
- The device should only be used for its intended purpose.
 - Please follow the procedures described in this instruction manual.
 - Keep the product away from children.
 - Do not use the device on damaged, irritated or diseased skin.

TIP

- The device must be checked for completeness of delivery and any visible damage.
- In the event of an incomplete delivery or damage due to faulty packaging or transport, please contact the service hotline.
- Before use, cleanse your facial skin and open your pores, e.g. with a warm towel or steam.

Technical data of the device

Symbol	24172
EAN	5907451335680
Material	plastic
Power supply	battery-powered
Landing	USB cable
Working time	1-2 hours
Charging time	2 hours
Device power	2.5W
Number of power levels	3

Startup/Use

- Please remove the packaging materials and dispose of them properly and clean the product and accessories.
- Before first use, charge the device using the USB cable. Connect the USB cable to a USB power source and connect the USB cable plug to the product's USB input port.
- Before the treatment, thoroughly cleanse your face and remove any makeup. For best results, it's recommended to open your pores, for example, with warm water or steam.
- Select the appropriate tip and attach it correctly to the device.
- Turn on the device by pressing the power button.
- Set the suction power to the lowest level, especially when using for the first time.
- Gently place the tip on your skin and move the device slowly across your face. Do not hold the device in one place for more than 2-3 seconds.
- Work from the center of the face outward, avoiding the eye and mouth area.
- If necessary, gradually increase the power level to suit your skin type.
- After completing the treatment, turn off the device.
- After completing the treatment, turn off the device.
- Remove the tip and clean it with a damp cloth or lukewarm water, then dry it thoroughly.

Safety tips and rules

- Do not use a product that may be damaged or has sharp edges or other dangerous areas that could cause injury.
- This product may be used by children aged 8 years and over and by persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge, provided they are supervised or have been instructed on the safe use of the product and understand the resulting risks.
- Always place the product on a stable and dry surface.
- For cleaning, use a damp cloth or a mild detergent.
- Do not use a damaged product.
- Do not use the product near heat sources, open flames or in damp rooms.
- The product is not intended for use by children.
- Cleaning and user maintenance must not be performed by children without supervision.
- The product itself and its accessories are not toys and should be kept out of the reach of children to avoid danger.
- Keep packaging and all protective materials out of the reach of children as they may pose a risk of suffocation.
- Protect the product from extreme temperatures, direct sunlight, strong vibrations, high humidity, moisture, flammable gases, vapors and solvents.
- Do not expose the product to mechanical stress.
- If safe operation is no longer possible, discontinue use and secure the product against further use. Safe operation is no longer possible if the product: - has been damaged, - is not functioning properly, - has been stored for an extended period under unfavorable conditions, or - has been subjected to excessive stress during transport.
- Do not use the product if any part is damaged. If the cord is damaged, do not attempt to repair it yourself.
- Do not hit or drop this product. Stop using it if it is damaged by a fall or other mechanical damage.
- The product should be stored in a cool and dry place.
- Make sure not to damage the USB cable with sharp edges or hot objects. Fully unwind the USB cable before use.
- The included USB cable is only suitable for this product.
- Turn off the product immediately and remove the cable from the product if it smells burning or smoke is visible. Have the product checked by a qualified technician before using it again.
- If you're using a USB power adapter, the power outlet must be easily accessible at all times so that you can easily unplug the USB power adapter in an emergency. Please also follow the manufacturer's operating instructions.
- Before each charging process, check the USB cable and device for any damage. Never use or charge the device if you notice any damage.
- The device is intended for external use on the skin only.
- Do not use the device on damaged, irritated, burned or infected skin.
- Do not use the device around the eyes, mouth or on very thin and sensitive skin.
- When using for the first time, always select the lowest suction power level.
- Do not hold the tip in one place for more than 2-3 seconds to avoid irritation or bruising.

- Do not use the device more than 2-3 times a week.
- If you experience pain, burning or discomfort, discontinue use immediately.
- The product should always be used as intended.

Safety instructions for batteries/rechargeable batteries

- Keep batteries/rechargeable batteries out of the reach of children. If swallowed, seek medical attention immediately!
- Disposable batteries must not be recharged. Batteries/rechargeable batteries should not be short-circuited and/or opened. This may result in overheating, fire, or explosion.
- Never throw batteries/rechargeable batteries into fire or water.
- Batteries/rechargeable batteries should never be subjected to mechanical stress.
- Risk of acid leakage from batteries/accumulators.
- Avoid extreme conditions and temperatures that may affect batteries/accumulators, e.g. radiators/direct sunlight.
- If batteries/rechargeable batteries leak, avoid contact with skin, eyes, and mucous membranes! Immediately flush the affected areas with clean water and seek medical attention!
- Leaking or damaged batteries/rechargeable batteries can cause chemical burns when in contact with skin. Therefore, wear appropriate protective gloves.
- The product has a built-in battery that cannot be replaced by the user.

Warning and information symbols

	<u>TIPS AND INFORMATION ON THE MANAGEMENT OF USED PACKAGING</u> The packaging is made of environmentally friendly materials that can be disposed of at your local recycling center. Used packaging material should be delivered to a designated waste disposal site designated by local authorities. For information on recycling options for used products, please contact your local municipality or city office. <u>DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT</u> For environmental reasons, used electrical and electronic products should not be disposed of with regular household waste, but should be disposed of properly. Information on collection points and their opening hours can be obtained from the relevant authorities. This product complies with the requirements of the applicable European and national directives. The product meets European and national requirements for the safety of devices and products. This product complies with the requirements of the relevant European and national RoHS directives.
---	---

	<u>Improper disposal of batteries/accumulators poses a threat to the environment!</u> Batteries/accumulators should not be disposed of with household waste. They may contain harmful heavy metals and should be treated as special waste. The chemical symbols for heavy metals are: Cd = cadmium, Hg = mercury, Pb = lead. Therefore, used batteries/accumulators should be taken to municipal hazardous waste collection points. The built-in battery must not be disassembled for disposal. The product should be taken to a collection point for used electronics.
---	--

We reserve the right to make changes to the text, design and product data without notice.

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